

SIEMENS

POLYMOBIL 10

SP

Planning Guide

© Siemens AG 1994

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Print No.: RXR8-120.021.01.03.02

Replaces: RXR8-120.021.01.02.02

English

Doc. Gen. Date: 04.99

Chapter	Page	Revision
all	all	03

1	General Notes	1 - 1
	General notes	1 - 1
	Safety	1 - 2
	System configuration	1 - 3
2	Room Planning	2 - 1
	Installation notes	2 - 1
	Dimensions: POLYMOBIL 10	2 - 2
3	Preparation for Installation	3 - 1
	Electrical installation on site	3 - 1
4	Technical Data	4 - 1
	Electrical data	4 - 1
	Weight and heat dissipation	4 - 1
	Environmental conditions	4 - 1
	Packing and transport routes	4 - 2
	Surface colors	4 - 2
5	Project Management	5 - 1
	Responsibility of the project manager towards the service contractor	5 - 1
6	Changes to Previous Version	6 - 1

This page intentionally left blank.

General notes



- With distribution of these revision level, all preceding planning guides, Speed - Infos (PG's) and drafts lose their validity.
- All layouts issued by the Planning Departments must bear a note referring to the installation and delivery conditions of Siemens Medical Engineering Group. The installation and delivery conditions must be submitted with the layouts.
- Unless otherwise specified, all dimensions are indicated in "mm".



- The symbol indicates a change (see revision status).



- Orientation points
Points specific to system components to which reference is made when positioning system components to each other or in the room.
The isocenter of a radiographic system is always illustrated as the orientation point.

- Fixpoints
Clearly marked points on system components, installation ceiling, walls or floor on which cable outlets are located.
Illustration in the drawings: circle with letter/number-combination.
The cable lengths establish the maximum fixpoint distances and thus the maximum distances between the individual system components.

- Room height
The room height is the distance measured from the top surface of the floor to the bottom surface of the ceiling structural elements (Unistrut rails) (bottom surface of drop ceiling).



- Room lighting
For the constancy test according to DIN 68 68-57, it must be possible to dim the lighting in rooms in which diagnosis is made on image playback devices (monitors).
For the constancy test, it is also required that the value documented in the acceptance test must be exactly reproducible.
- According to DIN 68 68-57 (international standard in preparation), the lighting in rooms in which image playback devices (monitors) are used for diagnosis, the following requirements must be met:
 - adjustable, no anti-glare screen, reproducible adjustment of the lighting (e. g. dimmer with scala),
 - no glare or reflection from windows, lights and light boxes in the standard working position of the monitors.


Hotline + 49 (9191) 18 - 8080

Safety

- The provisions of the relevant fire protection regulations must be observed for the premises.
- The system has been developed according to EN 60601 - 1.
- Minimum dimensions (e. g. room heights, safety distances) indicated in the planning guides are marked "min."
- ◆ - Basic strength against electromagnetic sources of interference.
Result of lightning discharges.
The protection targets of the different lightning protection areas up to the unit connection are also specified in the IEC 1024, DIN 48810, VDE 0675 and in the DEMVT recommendations.

System configuration

(schematic)

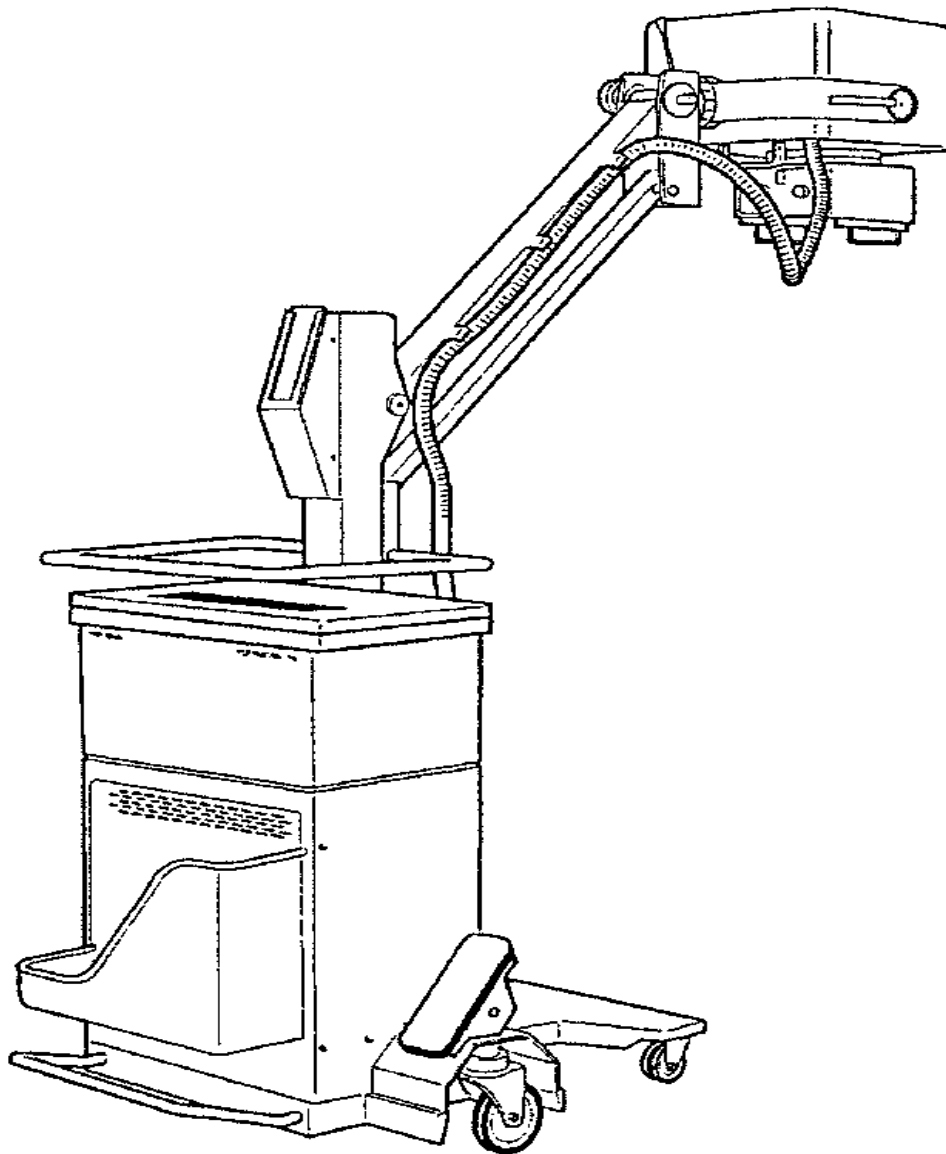


This page intentionally left blank.

Installation notes

- In the FRG (Germany), the POLYMOBIL 10 may be used in rooms intended for medical purposes only if the installation complies with DIN VDE 0107.
- In all countries outside the FRG (Germany), compliance with all country-specific regulations must be ensured.
- When planning the installation, ensure that the system's electrical plug can be disconnected from the mains quickly and without obstacles if a hazardous situation should arise.
- The mains cable (3 x 1.5²) is 5.5 m long and equipped with a 220 V plug.
- In countries where the 220 V plug is not permitted, it must be replaced with the plug required by the respective country.


Dimensions: POLYMOBIL 10

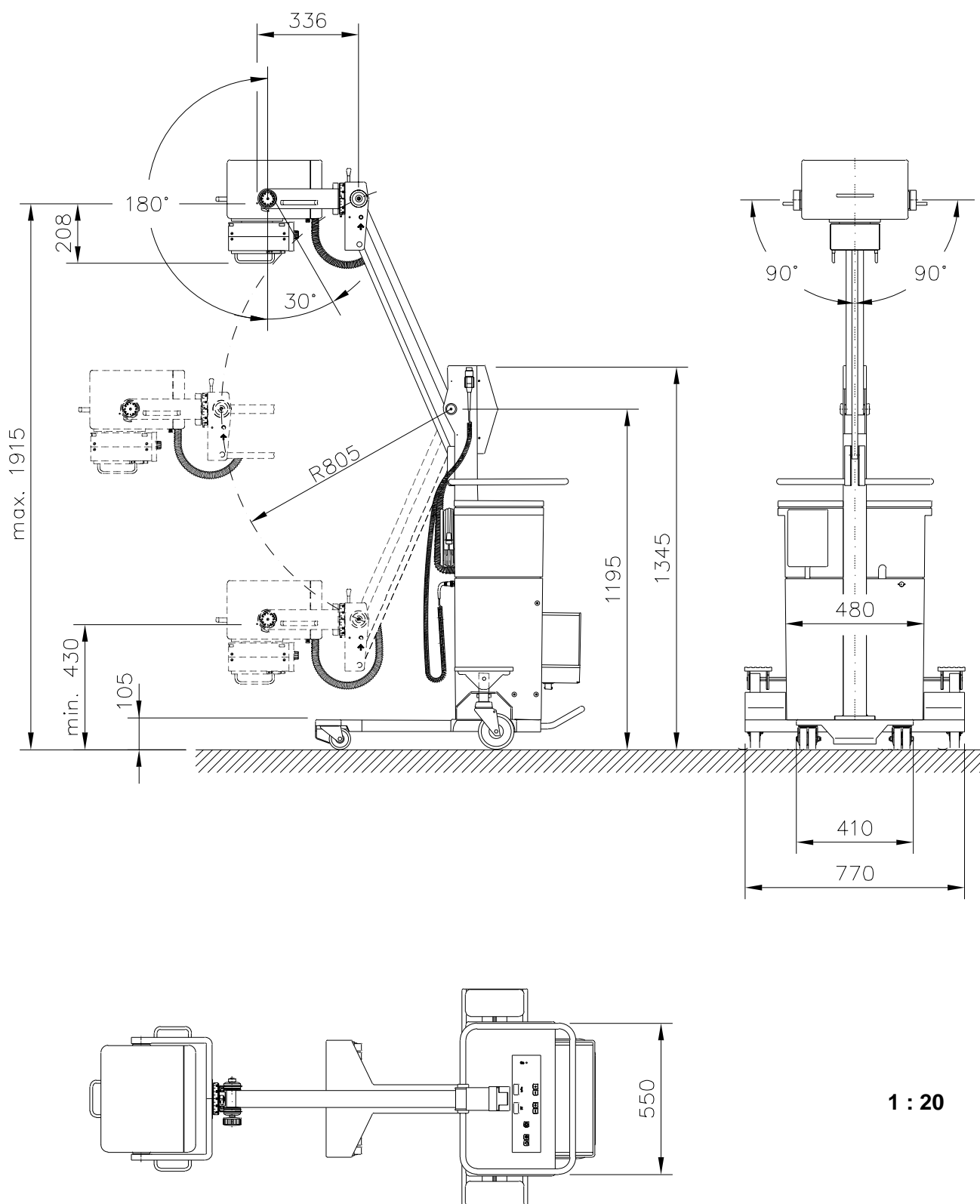

1 : 20

Electrical installation on site

Recommendation for on-site power distributor according to DIN VDE 0107

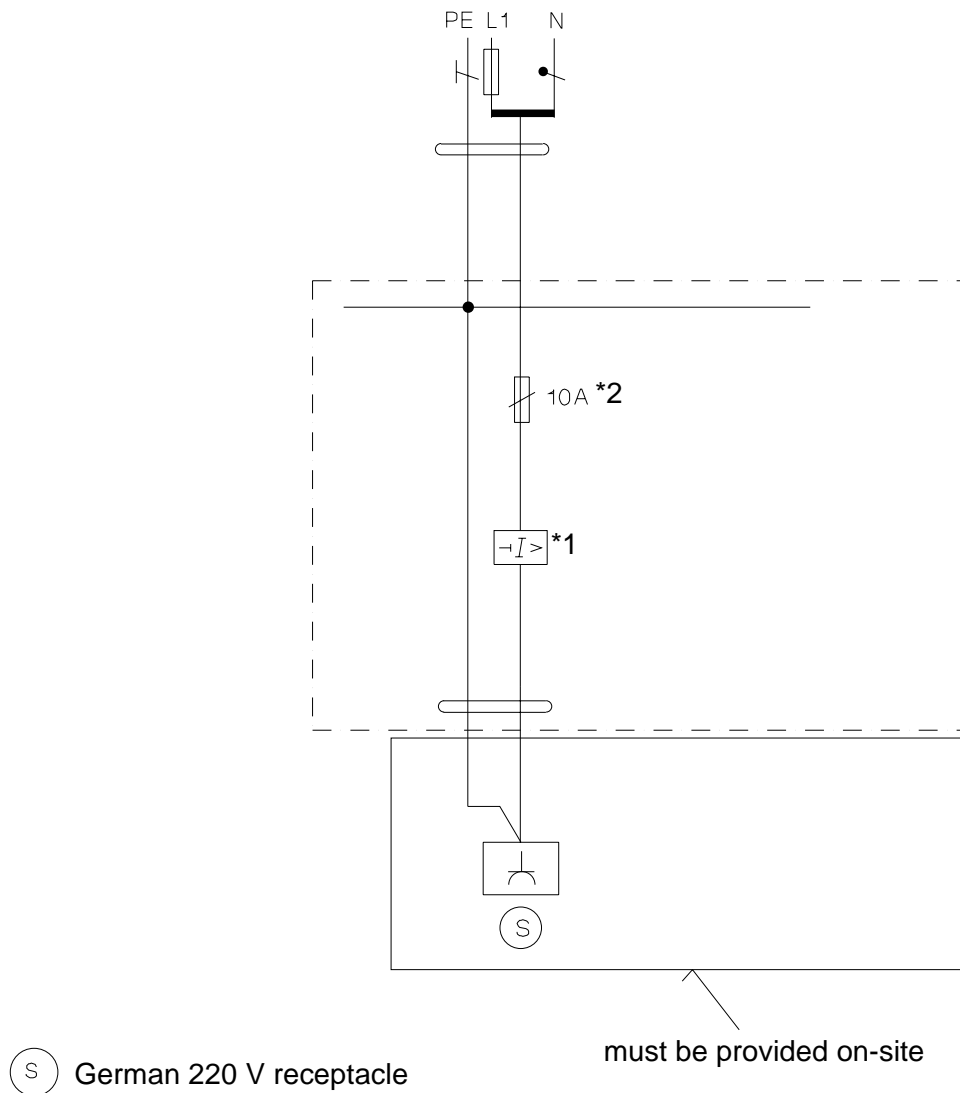

(S) German 220 V receptacle

*1 Fault current circuit breaker $I_{\Delta N}$ 30 mA for alternating and pulsating direct currents

*2 At 110 V / 16 A, at 230 V / 10 A

This page intentionally left blank.

Electrical data ♦

Mains connection	Internal line impedance	Fuse (internal)	Power consumption	Connection value according to DIN 6814 as the basis for calculating rates
1~100-143 V \pm 10 % 50/60 Hz 1~180-240 V \pm 10 % 50/60 Hz	max. 1.0 Ω	10 A	approx. 960 VA approx. 1776 VA	Formula: Mains voltage x fuse rating (internal)

Weight and heat dissipation

—	Weight [kg]	Heat dissipation in stand-by mode [W]
Total weight without packaging	approx. 173 kg	200 W
Total weight without packaging	approx. 245 kg	—

Environmental conditions ♦

POLYMOBIL 10	Operation	Transport / Storage
permissible ambient temperature	+ 10° C ... + 40° C	- 20° C ... + 70° C
permissible relative humidity	30 % ... 75 %	10 % ... 100 %
permissible air pressure	700 hPa ... 1060 hPa	500 hPa ... 1060 hPa

Packing and transport routes

Maximum size of crate	L 840 x W 1380 x H 1620 mm
Weight	approx. 245 kg
Door width required (without packaging)	800 mm

Surface colors

Main color	White pebbled, Med surface No. 4146 similar to RAL gray - white 9002
-------------------	---

Responsibility of the project manager towards the service contractor

The scope of the project manager's responsibilities requires that he

- is at the installation site when the system arrives
- supports the service contractor in solving problems
- clarifies the final location of the individual components
- checks to ensure that the installation is proceeding as specified
- clarifies problems together with the service contractor prior to the delivery of the system, e. g.
 - establishes the transport route of the truck
 - establishes the transport route within the building

NOTICE



The supervising SIEMENS project manager is responsible for the entire project management.

Furthermore, he is responsible for perfect and proper installation of the system.

Perform the further work according to the technical documentation (customer service instructions, installation instructions, etc.).

This page intentionally left blank.

Chapter	Page	Change
0-6		Layout change, thus Rev. - Stand of the document changed from 02 to 03
1	1-1 and 1-2	General Notes updated
4	4-1	Technical Data updated
5	5-1	Project Management updated
6	6-1	Changes to Previous Version updated

This page intentionally left blank.